

m.pen^[pro]
ultimate
microneedling
by mesoestetic®

85

Me

**perfect synergy
between technique
and drug substances**


Microneedling significantly increases the transepidermal penetration of drug substances¹

Can be combined with other medical-cosmetic techniques to deliver better results.


A very safe technique for all types of skin.

microneedling: What is it?


Microneedling is a minimally invasive technique that creates thousands of epidermal microcanals that allow substances to pass through the horny layer of the skin, while simultaneously producing dermoepidermal biostimulation, incrementing the activity of the drug substances.

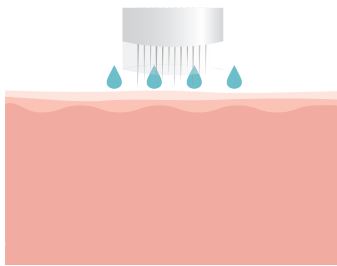

Many micropunctures are made, which in turn generate thousands of microcanals.

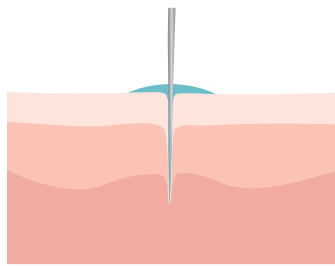

This facilitates the penetration of drug substances that are applied topically at the same time as microcanal creation.

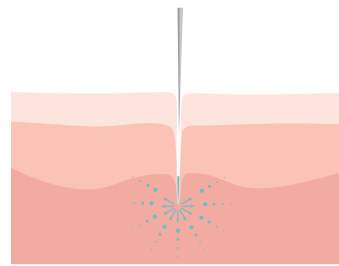

The drug substances reach the deepest skin layers, incrementing their effect.

Effects

01

Stimulates fibroblasts and the production of elastin and collagen I and III

02

Activates the skin repair mechanisms

03

Reduces expression lines and pore size

04

Reduces the appearance of hypertrophic scars and stretch marks

05

Improves skin firmness


06


Stimulates the hair follicles

m.pen^[pro] the ultimate microneedling solution


m.pen^[pro] is the new device developed to afford spectacular results safely and effectively. Headpiece movement is automatic, constant and perpendicular to the skin, contributing to its recovery.

What are its features?

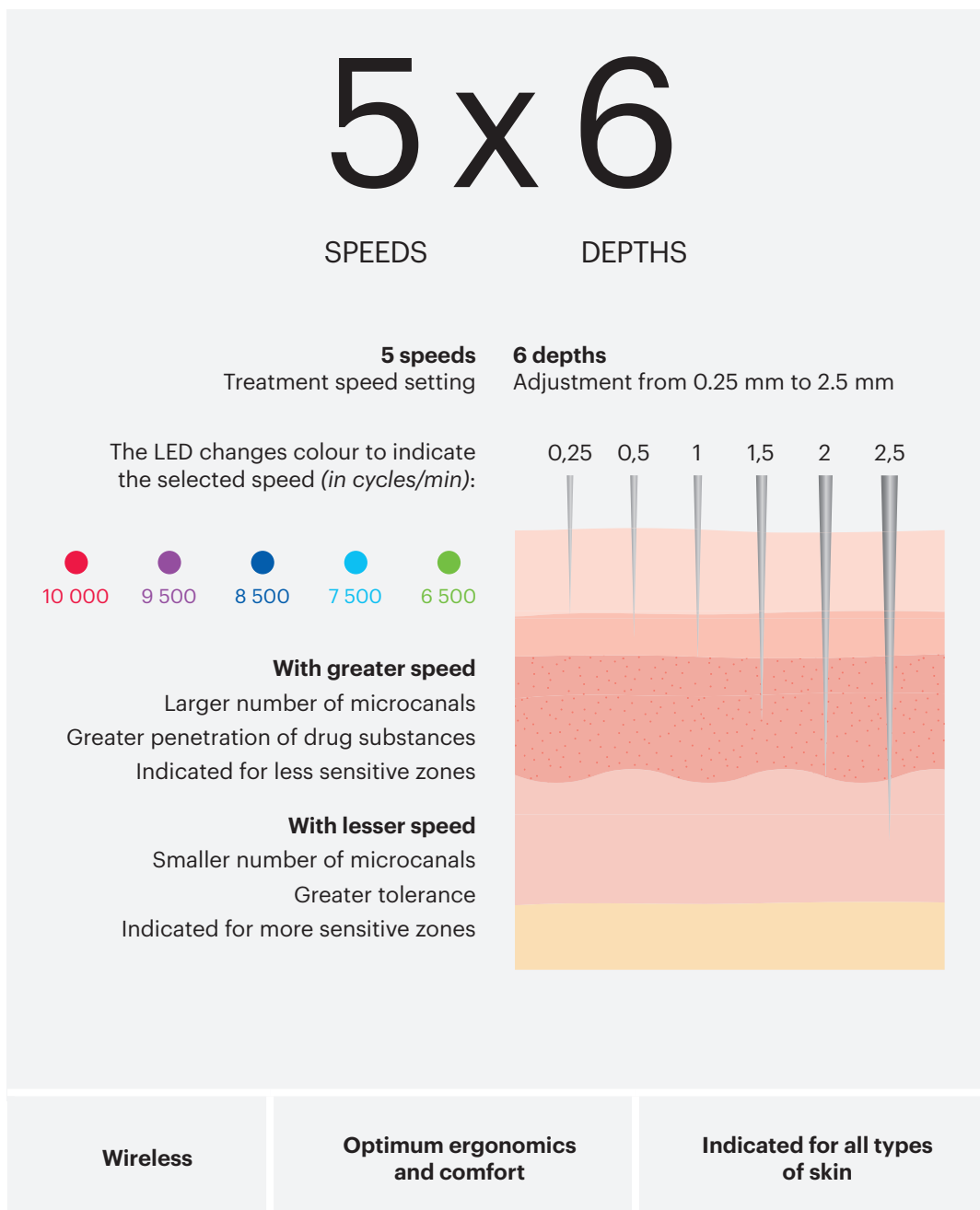

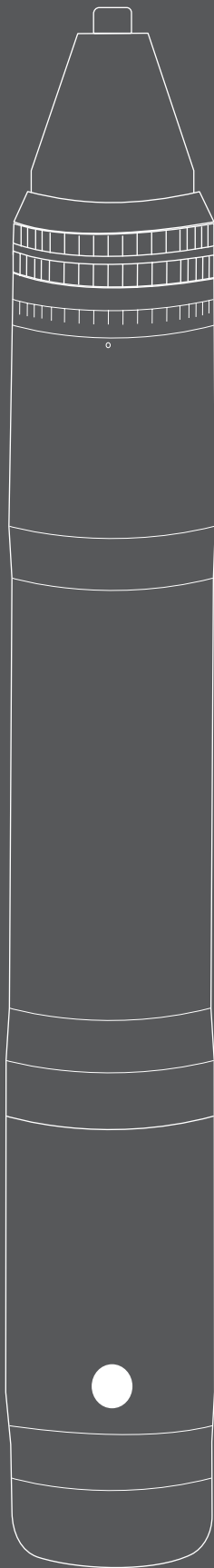


headpiece
zenithal
view

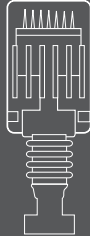

headpiece
lateral view

> Each headpiece consists of 12 microneedles with a circular configuration, allowing improved access to small zones and folds.





> Single-use sterile headpieces supplied in packs of 10 units, each presented individually to avoid the risk of contamination.

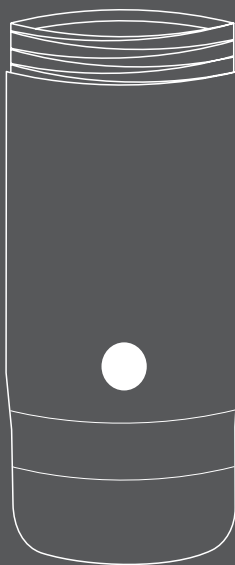

> Two exchangeable batteries that each provide 2 hours of use, enabling uninterrupted treatment and simultaneous charging during the procedure.

> The batteries can be charged by plugging the adaptor into the mains or through a USB port.

Application zones and regimens

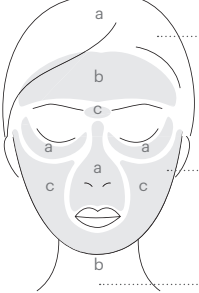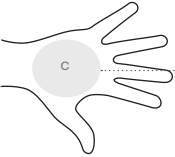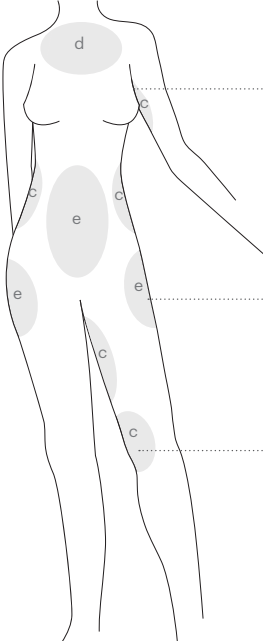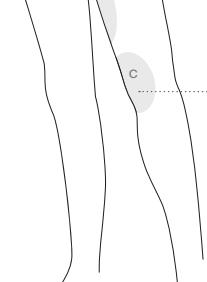APPLICATION ZONES*	INDICATIONS
 scalp face, neckline, hands and other body zones neck, face and double chin	hair loss hyperpigmentations wrinkles
 hands, face and neckline	lack of vitality
 inner surfaces of the arms and legs, abdomen and flanks hip fat, double chin, flanks and abdominal fat	flaccidity and loss of turgidity localised fat
 knees (inner surface) and legs	cellulite
* Application zones and relation to microneedle depth: a. 0.25 mm b. < 0.50 mm c. < 1 mm d. < 1.5 mm e. > 1.5 mm	

m.pen [pro] can be displaced along a straight line or transversely, and allows the tracing of small circles. The technique varies according to the treatment zone and programme.

The efficacy of the synergy between **m.pen [pro]** and **c.prof** and **x.prof** has been demonstrated clinically under medical supervision.

RECOMMENDED PRODUCT	SESSIONS	TIME BETWEEN SESSIONS
c.prof 230 hair loss solution	5	15 days
c.prof 210 depigmentation solution	5	15 days
c.prof 213 mesotox solution	5	15 days
c.prof 211 photoaging solution	5	15 days
c.prof 222 body firming solution	10	7 days
c.prof 221 lipolytic solution	10	7 days
c.prof 220 draining solution	10	7 days



Protocol

Pre-treatment recommendation

It is advisable to apply a topical anaesthetic onto the treatment zone 30-60 minutes before the session.

01/ Cleaning

Clean the skin with **hydramilk** and **hydratonic**, and disinfect the zone.

hydramilk

500 ml

hydratonic

500 ml



02/ Application of the product and technique

Apply the solution chosen from the **c.prof** or **x.prof** product range onto the treatment zone, combined with the use of the **m.pen** ^[pro]

c.prof/x.prof

m.pen ^[pro]



03/ Application of the technique

Perform the technique on the skin following the tracings in different directions.

m.pen^[pro]



04/ Soothing and repairing

After completing the session, apply the **crystal fiber mask** to soothe, repair and decongest the skin.

crystal fiber mask
5 x 25 ml



Post-treatment recommendation

**moisturising
sun protection**
500 ml




Complete the treatment by applying **moisturising sun protection SPF 50+** to protect the skin against the harmful effects of UV radiation.

The sessions are performed once a week or once every two weeks, depending on the treatment zone and the treatment aims, and always according to professional judgement.

Each session lasts between 30 and 60 minutes, and the patient can resume their daily living immediately.

Demonstrated efficacy


In vivo histological view showing the microcanals formed with the microneedling technique¹.

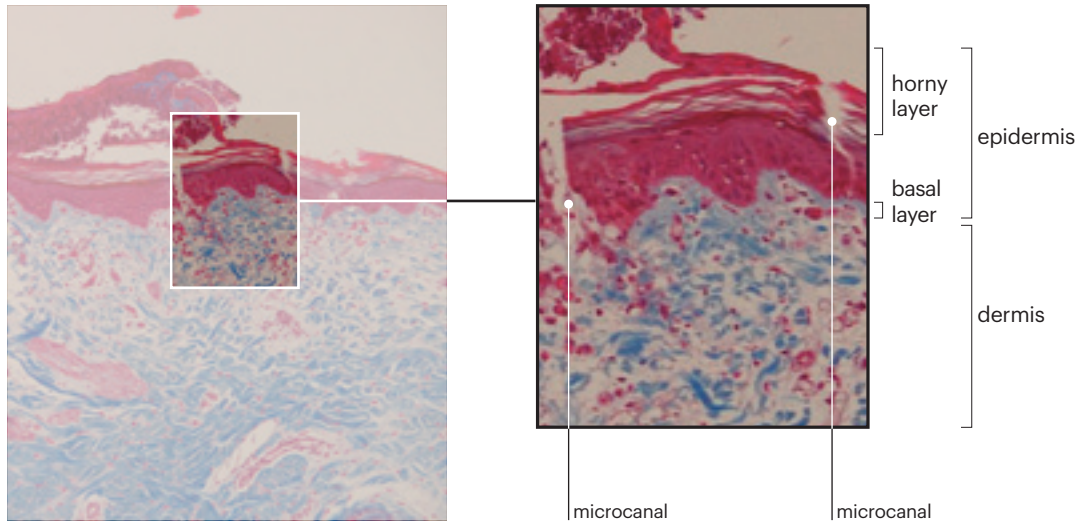


Stretch marks

m.pen^[pro] + **x.prof 040** centella asiatica, **x.prof 014** 1% organic silicon, **x.prof 022** piruvex



sex: female **age:** 49 years **SPEED** > 4th level (9500 cycles/min) **DEPTH** > 1.5 mm

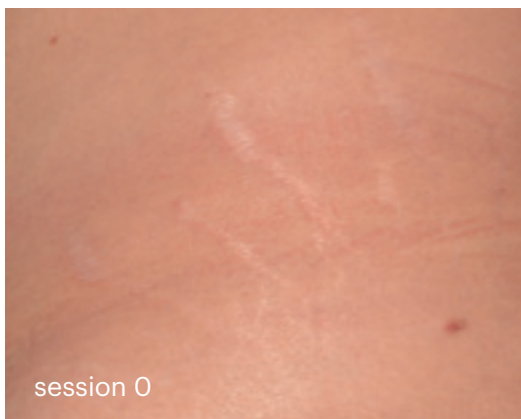

A. stretch mark size 6 mm
B. dermis thickness 1.2 mm

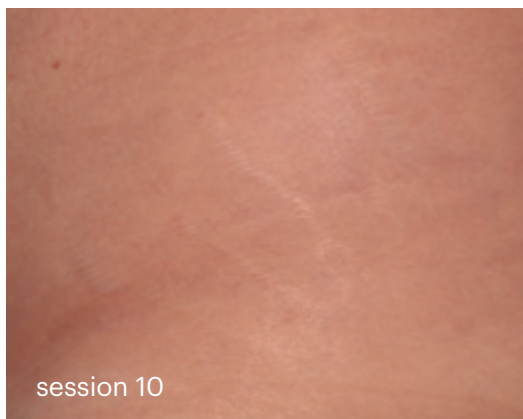

A. stretch mark size 3.4 mm (- 43%)
B. dermis thickness 1.7 mm (+ 41%)

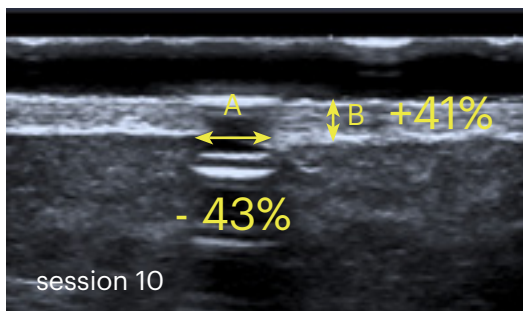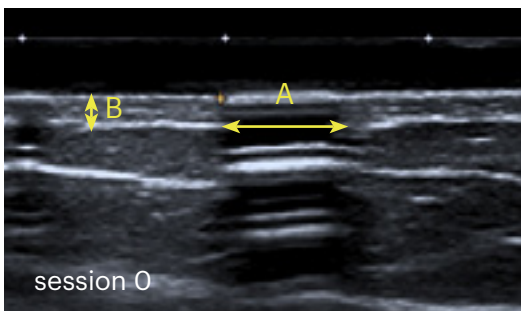

Source: Department of Medical Affairs, Mesoestetic Pharma Group

Double chin

m.pen^[pro] + **x.prof 014** 1% organic silicon, **x.prof 012** l-carnitine and **x.prof 015** melilot and rutin extract

sex: female **age:** 47 years **SPEED** according to tolerance **DEPTH** < 1 mm




Source: Department of Medical Affairs, Mesoesthetic Pharma Group



Depigmentation

m.pen^[pro] + **c.prof 210** depigmentation solution

sex: female **age:** 33 years **SPEED** > 3rd level (8500 cycles/min) **DEPTH** > 1.5 mm

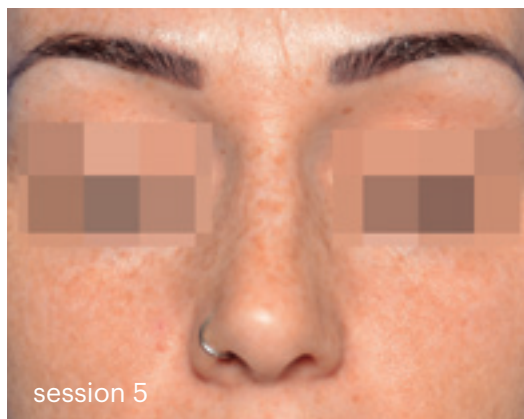

Appearance under UV light

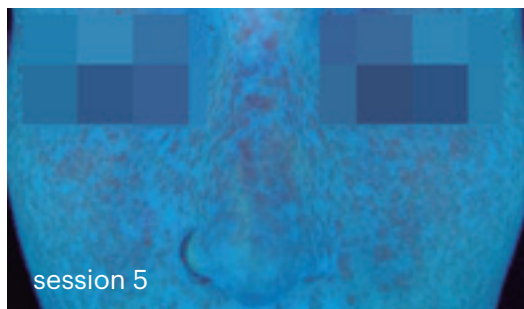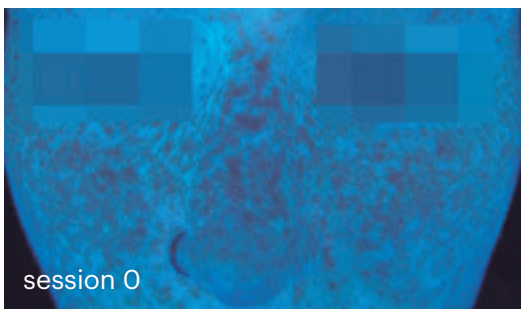

Source: Department of Medical Affairs, Mesoesthetic Pharma Group

m.pen^[pro]

The definitive advance in the field of dermoepidermal biostimulation.

Effective, rapid and safe transepidermal device.

5 x 6. Adjustment of micropuncture depth and speed for maximum treatment personalisation.

Reinforces treatment efficacy.

Clinically demonstrated efficacy.

Visible results in only 3 sessions.



The pack includes:

m.pen^[pro] device, 2 headpieces, 2 exchangeable batteries, adaptor + power cable.
The headpieces can be acquired individually in packs of 10 units each.

1. Ignacio Ordiz: "Transepidermal drug administration routes". Doctoral Thesis, University of Oviedo, 2016. <https://www.educacion.gob.es/teseo/mostrarRef.do?ref=1276689>



mesoestetic Pharma Group, s.l.
C/ Tecnologia, 25
08840 - Viladecans (Barcelona)
Tel. +34 93 325 20 30
www.mesoestetic.com
made in Spain

KFEPN0003/2212